




Weck™ Hem-o-lok™ Purple*Plus*™

Endoscopic Polymer Clip Ligation System

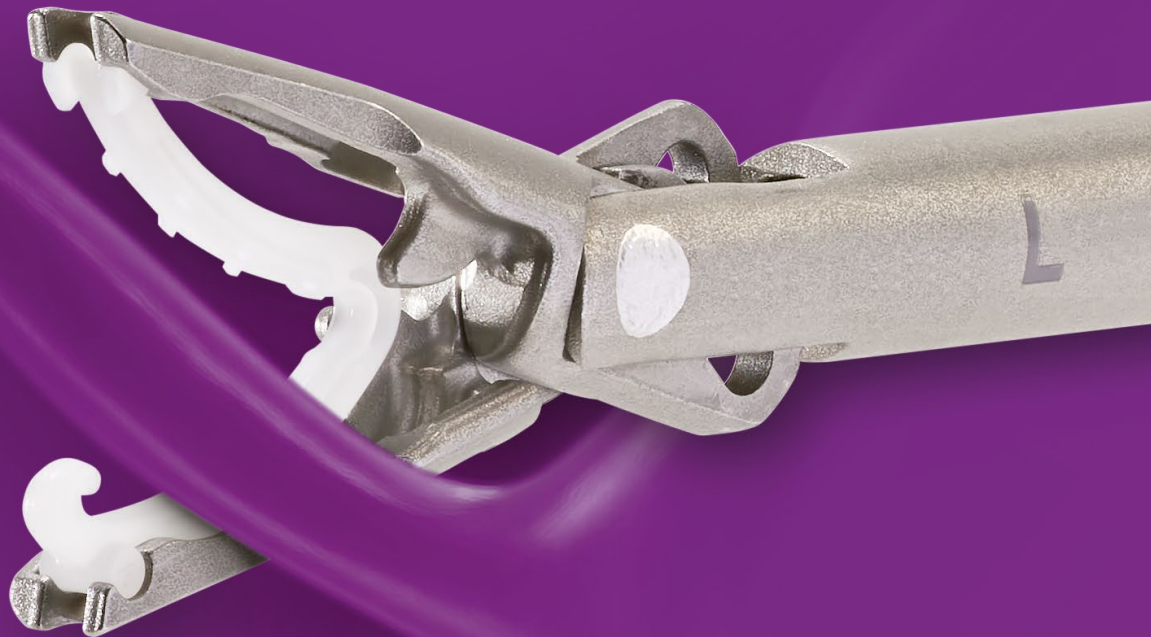

NEXT GENERATION SYSTEM

Hem-o-lok™ Purple*Plus*™ Polymer Clip Ligation System now offers enhanced clip stability and procedural flexibility with the trusted security of Weck™ Hem-o-lok™ Clips.

Teleflex™
Empowering the future of healthcare

Less space. More possibilities.

Redesigned for **improved clip stability** during clip application and the flexibility to ligate 2–13 mm vessels through a **small 5 mm trocar**.

Ligate larger vessels through a 5 mm port


The flexible hinge of Hem-o-lok™ PurplePlus™ Large Polymer Clips allow ligation of 2–13 mm vessels and tissue structures through a 5 mm port, further advancing the minimally invasive approach without compromising jaw aperture.

Designed to prevent clip fall-out

Hem-o-lok™ PurplePlus™ Large Polymer Clips offer three points of connection within the applier jaws, so clips remain secure from loading to latching on tissue structures and vessels.

Secure closure you can feel and hear

Hem-o-lok™ PurplePlus™ Large Polymer Clips feature distinct tactile and auditory feedback to provide assurance of secure closure on targeted tissue structures and vessels.

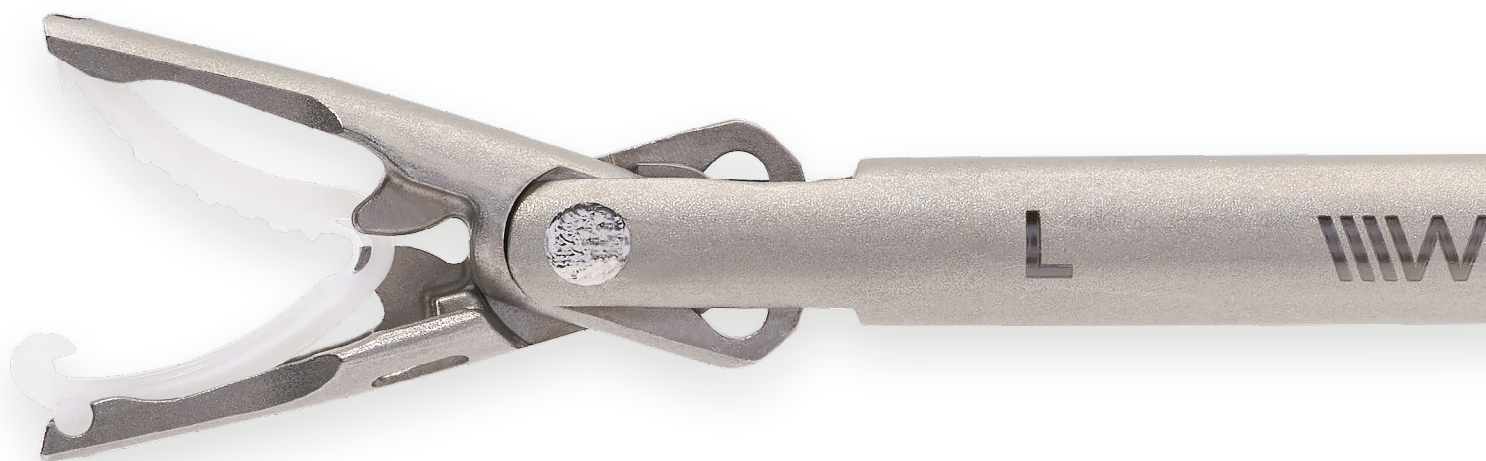

weck™
Hem-o-lok™
Endoscopic Polymer Clip Ligation System


PurplePlus™



eck®

5MM



CLINICIAN

Designed to give surgeons peace of mind, its 3 points of contact in the applier jaws **ensures a stable application on tissues and vessels.**



INSTITUTION

Hem-o-lok™ Purple**Plus**™ Polymer Clips are designed with 3 points of connection in the applier jaws, **reducing clip fall-out, minimizing waste, and maintaining patient safety.**



PATIENT

Ligation of larger tissue structures through smaller ports **helps to minimize patient trauma.**

Designed for confidence

Hem-o-lok™ PurplePlus™ Large Polymer Clips are designed to provide secure closure for peace of mind during surgical procedures that involve larger vessels and tissue structures.



Hem-o-lok™ PurplePlus™ Clip Cartridge features


- Designed for **easy and precise loading** that can enhance surgical procedure efficiency
- **Arrow and hook icons** to clearly indicate clip orientation
- **Upgraded proprietary branding, including yellow accents**, to ensure you are using genuine Weck™ Cartridges and Clips



PurplePlus™

Hem-o-lok™ PurplePlus™ Large Polymer Clip features

- With **3 points of contact** in the applier, Hem-o-lok™ PurplePlus™ Large Polymer Clips remain secure during loading and latching, resulting in a lower chance of dropping clips during loading and latching.
- The **bosses** and third point of contact of the Hem-o-lok™ PurplePlus™ Polymer Ligating Clip enhance security in the applier jaws.
- **Bow-shaped design** of the clip allows for easy removal with the corresponding Hem-o-lok™ PurplePlus™ Endoscopic Remover.
- The Hem-o-lok™ PurplePlus™ Polymer Ligating Clip features **teeth and a locking design** that prevents slippage, providing a reliable 360-degree grip around a vessel or tissue structure.
- With its **flexible hinge**, Hem-o-lok™ PurplePlus™ Large Polymer Clips can be inserted through 5 mm ports to ligate 2–13 mm vessels and tissues, reducing the need for additional and larger ports.
- The **locking mechanism** of Hem-o-lok™ PurplePlus™ Polymer Ligating Clip provides both tactile and auditory feedback, helping to eliminate second-guessing and doubts about secure closure.
- Made with **non-absorbable polymer**, Hem-o-lok™ PurplePlus™ Clips don't interfere with CT, MRI or X-ray diagnostics.

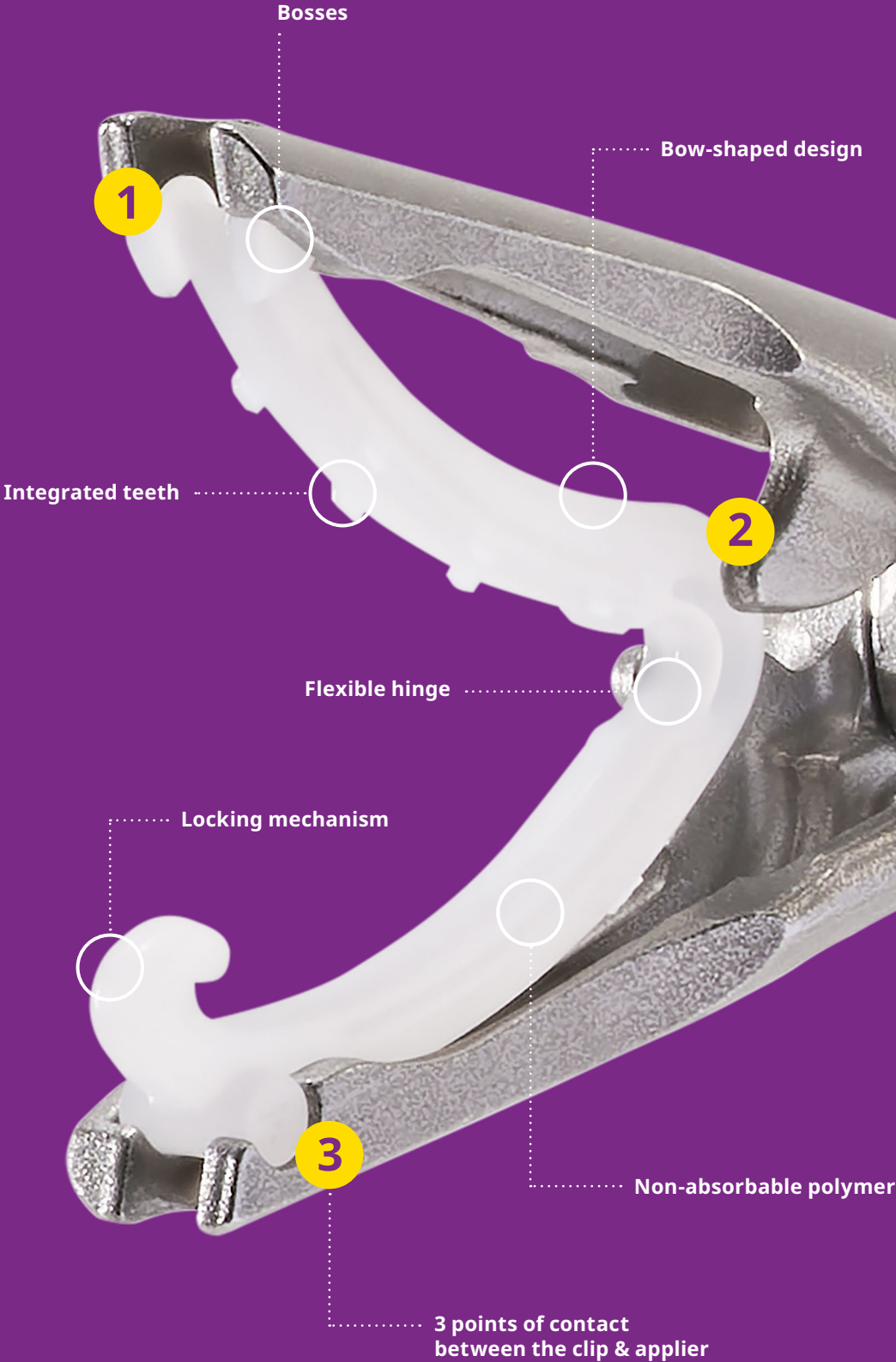

Unparalleled security*

Hem-o-lok™ Purple**Plus**™ Large Polymer Clips are engineered for unparalleled security with greater lateral clip retention and greater leak resistance during manual cold ligation on vessels and tissue structures than leading metal clip competitors.*



During lateral pull-off tests and pressure tests, Hem-o-lok™ Purple**Plus**™ Large Polymer Clips provided:

1.5X


Greater leak resistance*

14X

Greater clip retention*

Reinventing cold ligation, again

Our Hem-o-lok™ Purple**Plus**™ Polymer Clip Ligation System delivers a **large clip** with a smaller applicator—enabling an access profile that is **4X smaller** than our legacy 10 mm platform.

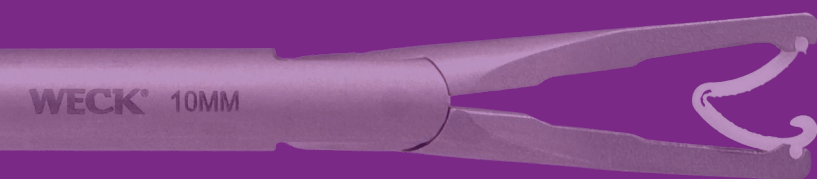

Legacy Hem-o-lok™ Large System / 10 mm platform



Hem-o-lok™ Purple**Plus**™ Large System / 5 mm platform



Hem-o-lok™ PurplePlus™ Polymer Ligating System (Ordering Information)

PRODUCT CODE	PRODUCT DESCRIPTION	QTY
HOL2LG	Hem-o-lok™ PurplePlus™ Large, Clips	6/CA 84/BOX
ME05LG	Hem-o-lok™ PurplePlus™ Large, Manual Endo, 5mm, 35cm, Clip Applier	1 Each
ME05RVLG	Hem-o-lok™ PurplePlus™ Large, Manual Endo, 5mm, 35cm, Clip Remover	1 Each

* Data on file at Teleflex, D054562. Test data may not be indicative of clinical performance. Compared to: Ethicon (Ligamax, LigaClip), Medline/Covidien (EndoClip)

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Rx Only

Contraindications: Weck™ Hem-o-lok **Plus**™ Ligating Clips are contraindicated for use as a fallopian contraceptive tubal occlusion device. Weck™ Hem-o-lok **Plus**™ Clips are contraindicated for use in ligating the renal artery during laparoscopic donor nephrectomies. The use of Weck™ Hem-o-lok **Plus**™ Ligating Clips with any other ligation system is contraindicated due to their incompatibility. The use of Weck™ Hem-o-lok **Plus**™ Endoscopic Ligating Clip Appliers and Removers with any other polymer ligating clips is contraindicated due to their incompatibility.

Teleflex, the Teleflex logo, Hem-o-lok, Hem-o-lok Plus, PurplePlus, and Weck are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Other names may be trademarks of their respective owners.

Information in this material is not a substitute for the product Instructions for Use. Refer to package insert for complete warnings, indications, contraindications, precautions, potential complications and instructions for use. Not all products may be available in all countries. Please contact your local representative. Revised: 02/2025

© 2025 Teleflex Incorporated. All rights reserved.

MC-011082 · REV 0